Leader in Biomolecular Solutions for Life Science

Phospho-NF-kB p65/RelA-S536 Rabbit pAb



Catalog No.: AP0475

57 Publications

Basic Information

Observed MW

65kDa

Calculated MW

60kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

NF-kappa-B is a ubiquitous transcription factor involved in several biological processes. It is held in the cytoplasm in an inactive state by specific inhibitors. Upon degradation of the inhibitor, NF-kappa-B moves to the nucleus and activates transcription of specific genes. NF-kappa-B is composed of NFKB1 or NFKB2 bound to either REL, RELA, or RELB. The most abundant form of NF-kappa-B is NFKB1 complexed with the product of this gene, RELA. Four transcript variants encoding different isoforms have been found for this gene.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

Immunogen Information

Gene ID5970

Swiss Prot
Q04206

Immunogen

A synthetic phosphorylated peptide around S536 of human NF-kB p65/RelA (NP_068810.3).

Synonyms

p65; CMCU; NFKB3; AIF3BL3; Phospho-NF-kB p65/RelA-S536

Contact

www.abclonal.com


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

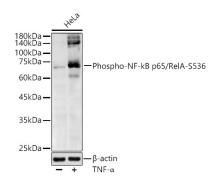

GAPDH TNF-α Western blot analysis of lysates from HeLa cells using Phospho-NF-kB p65/RelA-S536 Rabbit pAb (AP0475) at 1:1000 dilution. HeLa cells were treated by TNF- α (20 ng/mL) at 37°C for 30 minutes.


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.


Detection: ECL Basic Kit (RM00020).

Exposure time: 60s.

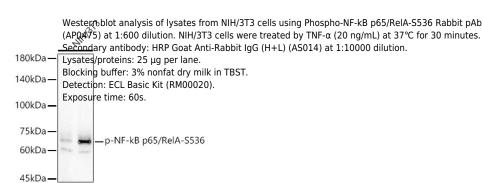

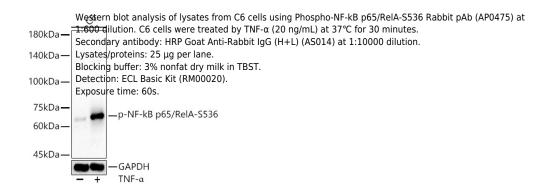

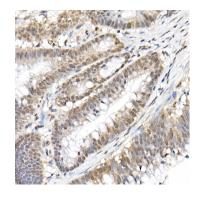

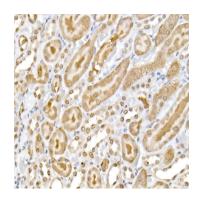

Validation Data

Immunohistochemistry analysis of Phospho-NF-kB p65/RelA-S536 in paraffin-embedded human colon carcinoma using Phospho-NF-kB p65/RelA-S536 Rabbit pAb (AP0475) at dilution of 1:50 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Immunohistochemistry analysis of Phospho-NF-kB p65/RelA-S536 in paraffin-embedded rat kidney using Phospho-NF-kB p65/RelA-S536 Rabbit pAb (AP0475) at dilution of 1:50 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.