

A11223

Leader in Biomolecular Solutions for Life Science



## Ku70 Rabbit mAb

Catalog No.: A11223

Recombinant

1 Publications

### Basic Information

#### Observed MW

70kDa

#### Calculated MW

70kDa

#### Category

SMab Recombinant Monoclonal  
Antibody

#### Applications

WB,IHC-P,IF/ICC,ELISA

#### Cross-Reactivity

Human,Mouse

#### CloneNo number

ARC0551

### Background

The p70/p80 autoantigen is a nuclear complex consisting of two subunits with molecular masses of approximately 70 and 80 kDa. The complex functions as a single-stranded DNA-dependent ATP-dependent helicase. The complex may be involved in the repair of nonhomologous DNA ends such as that required for double-strand break repair, transposition, and V(D)J recombination. High levels of autoantibodies to p70 and p80 have been found in some patients with systemic lupus erythematosus.

### Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

IF/ICC 1:50 - 1:200

### Immunogen Information

#### Gene ID

2547

#### Swiss Prot

P12956

#### Immunogen

A synthetic peptide corresponding to a sequence within amino acids 500-609 of human Ku70 (P12956).

#### Synonyms

ML8; KU70; TLAA; CTC75; CTCBF; G22P1; Ku70

### Contact



[www.abclonal.com](http://www.abclonal.com)

### Product Information

#### Source

Rabbit

#### Isotype

IgG

#### Purification


Affinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,0.05% BSA,50% glycerol,pH7.3.

## Validation Data

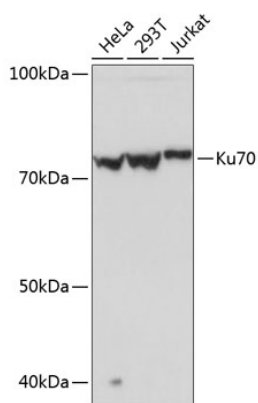

Western blot analysis of various lysates using Ku70 Rabbit mAb (A11223) at 1 : 1000 dilution.


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 1s.

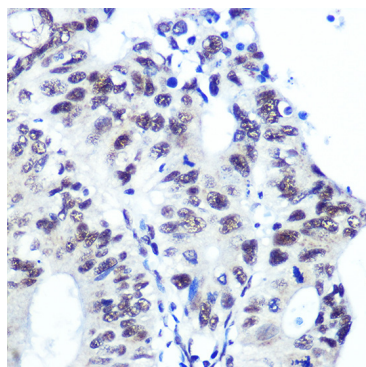

Immunohistochemistry analysis of paraffin-embedded human colon carcinoma using Ku70 Rabbit mAb (A11223) at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

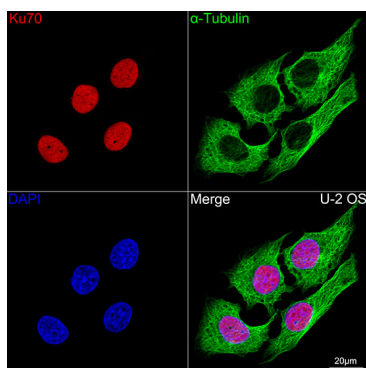

Confocal imaging of U-2 OS cells using Ku70 Rabbit mAb (A11223, dilution 1:100) (Red). The cells were counterstained with α-Tubulin Mouse mAb (AC012, dilution 1:400) (Green). DAPI was used for nuclear staining (blue). Objective: 100x.