A17991

Leader in Biomolecular Solutions for Life Science

MT-ND6 Rabbit pAb

Catalog No.: A17991 7 Publications



Basic Information

Observed MW 24kDa

Calculated MW 19kDa

Category Polyclonal Antibody

Applications WB,IHC-P,ELISA

Cross-Reactivity Human,Mouse,Rat

Background

Enables NADH dehydrogenase (ubiquinone) activity. Involved in mitochondrial electron transport, NADH to ubiquinone and mitochondrial respiratory chain complex I assembly. Predicted to be located in mitochondrial inner membrane. Implicated in Leber hereditary optic neuropathy; Leigh disease; and spinal muscular atrophy with lower extremity predominante 2B.

Recommended Dilutions

Immunogen Information

| 1:1000 - 1:5000 | Gene ID | Swiss Prot |
|-----------------|---------|------------|
| | 4541 | P03923 |
| 1:50 - 1:200 | | |

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 80-174 of human MT-ND6 (YP_003024037.1).

Synonyms MTND6; ND6; MT-ND6

€

WB

IHC-P

Product Information

www.abclonal.com

Isotype IgG Purification Affinity purification


Storage


Source


Rabbit

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.


Validation Data

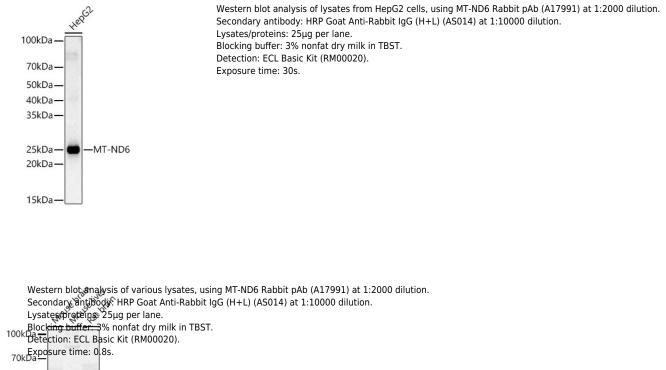

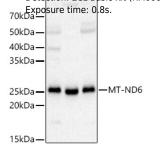


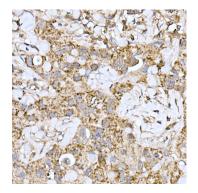

Immunohistochemistry analysis of MT-ND6 in paraffin-embedded human breast cancer using MT-ND6 Rabbit pAb (A17991) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

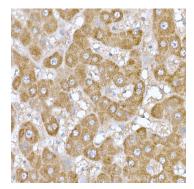

Immunohistochemistry analysis of MT-ND6 in paraffin-embedded human liver using MT-ND6 Rabbit pAb (A17991) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

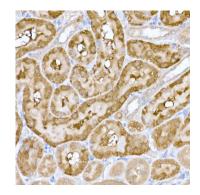

Immunohistochemistry analysis of MT-ND6 in paraffin-embedded rat kidney using MT-ND6 Rabbit pAb (A17991) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.