

A19710

Leader in Biomolecular Solutions for Life Science



SOX9 Rabbit mAb

Catalog No.: A19710

Recombinant

28 Publications

Basic Information

Observed MW

56kDa

Calculated MW

56kDa

Category

SMab Recombinant Monoclonal Antibody

Applications

WB,IHC-P,IF/ICC,IP,ELISA

Cross-Reactivity

Human,Mouse,Rat

CloneNo number

ARC0190

Background

The protein encoded by this gene recognizes the sequence CCTTGAG along with other members of the HMG-box class DNA-binding proteins. It acts during chondrocyte differentiation and, with steroidogenic factor 1, regulates transcription of the anti-Muellerian hormone (AMH) gene. Deficiencies lead to the skeletal malformation syndrome campomelic dysplasia, frequently with sex reversal.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:100 - 1:500

IF/ICC 1:50 - 1:200

IP 0.5µg-4µg antibody for
200µg-400µg extracts of
whole cells

Immunogen Information

Gene ID

6662

Swiss Prot

P48436

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 169-300 of human SOX9 (P48436).

Synonyms

CMD1; SRA1; CMPD1; SRXX2; SRXY10; SOX9

Contact



www.abclonal.com

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,0.05% BSA,50% glycerol,pH7.3.


Validation Data

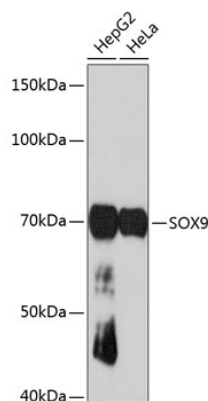


Western blot analysis of various lysates using SOX9 Rabbit mAb (A19710) at 1:1000 dilution.
Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25µg per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 3min.

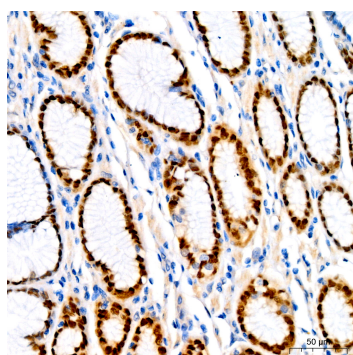


Immunohistochemistry analysis of SOX9 in paraffin-embedded human colon using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens). Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

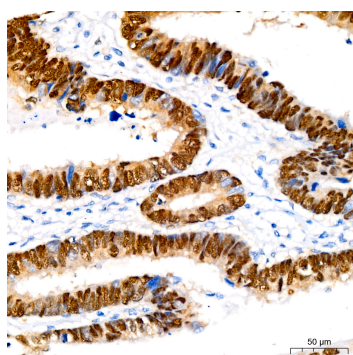


Immunohistochemistry analysis of SOX9 in paraffin-embedded human gastric cancer using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens). Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

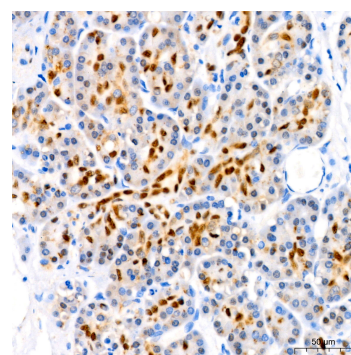


Immunohistochemistry analysis of SOX9 in paraffin-embedded human pancreas using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens). Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

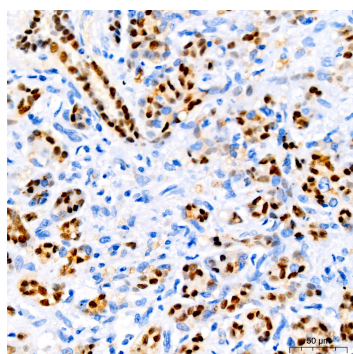

Immunohistochemistry analysis of SOX9 in paraffin-embedded Human pancreatic ductal carcinoma using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens). Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

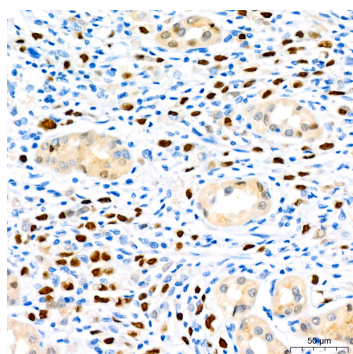


Immunohistochemistry analysis of SOX9 in paraffin-embedded human stomach using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens). Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

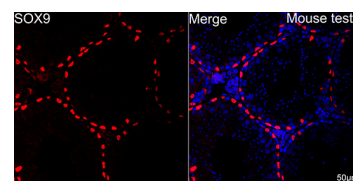


Confocal imaging of paraffin-embedded Mouse testis using SOX9 Rabbit mAb (A19710, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform high pressure antigen retrieval with 0.01M citrate buffer (pH 6.0) prior to IF staining.


Validation Data

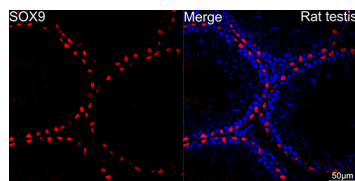

Confocal imaging of paraffin-embedded Rat testis using SOX9 Rabbit mAb (A19710, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform high pressure antigen retrieval with 0.01M citrate buffer (pH 6.0) prior to IF staining.

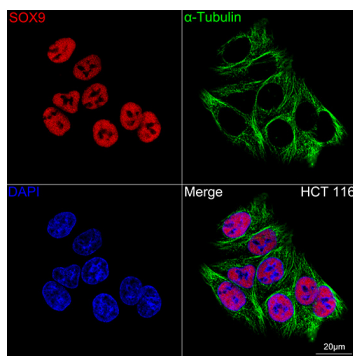

Confocal imaging of HCT 116 cells using SOX9 Rabbit mAb (A19710, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α-Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

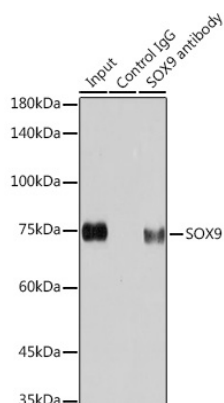

Immunoprecipitation analysis of 200 µg extracts of HeLa cells using 3 µg SOX9 antibody (A19710). Western blot was performed from the immunoprecipitate using SOX9 antibody (A19710) at a dilution of 1:1000.