## A19798

Leader in Biomolecular Solutions for Life Science

# VPS11 Rabbit mAb

Catalog No.: A19798 Recombinant



# **Basic Information**

Observed MW 108kDa

Calculated MW 108kDa

Category SMab Recombinant Monoclonal Antibody

Applications WB,IHC-P,ELISA

Cross-Reactivity Human, Mouse

CloneNo number ARC2325

# Background

Vesicle mediated protein sorting plays an important role in segregation of intracellular molecules into distinct organelles. Genetic studies in yeast have identified more than 40 vacuolar protein sorting (VPS) genes involved in vesicle transport to vacuoles. This gene encodes the human homolog of yeast class C Vps11 protein. The mammalian class C Vps proteins are predominantly associated with late endosomes/lysosomes, and like their yeast counterparts, may mediate vesicle trafficking steps in the endosome/lysosome pathway. Alternative splicing results in multiple transcript variants.

# **Recommended Dilutions**

### **Immunogen Information**

WB	1:500 - 1:1000	Gene ID
	1:50 - 1:200	55823
IHC-P	1:50 - 1:200	

Swiss Prot Q9H270

#### Immunogen

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human VPS11 (Q9H270).

#### Synonyms

END1; PEP5; HLD12; RNF108; hVPS11; DYT32; VPS11

Contact
---------

€

# **Product Information**

www.abclonal.com

**lsotype** IgG Purification Affinity purification


#### Storage


Source


Rabbit

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% sodium azide,0.05% BSA,50% glycerol,pH7.3.


# **Validation Data**

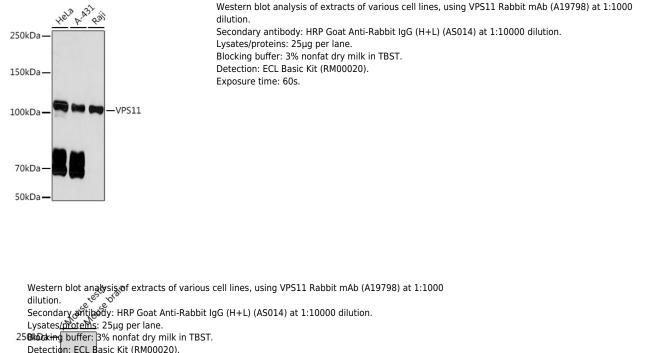

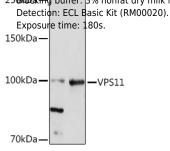

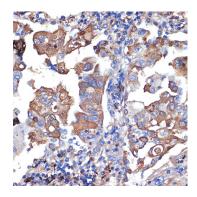

Immunohistochemistry analysis of paraffinembedded human lung cancer using VPS11 Rabbit mAb (A19798) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

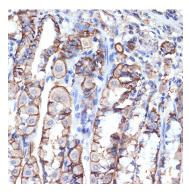

Immunohistochemistry analysis of paraffinembedded mouse stomach using VPS11 Rabbit mAb (A19798) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.