# ABclonal www.abclonal.com

# **GNAI2** Rabbit mAb

Catalog No.: A20930 Recombinant

## **Basic Information**

## **Observed MW**

40kDa

## **Calculated MW**

40kDa

#### Category

Primary antibody

### **Applications**

ELISA, WB, IHC-P

#### **Cross-Reactivity**

Human, Mouse, Rat

#### CloneNo number

ARC2911

# Background

The protein encoded by this gene is an alpha subunit of guanine nucleotide binding proteins (G proteins). The encoded protein contains the guanine nucleotide binding site and is involved in the hormonal regulation of adenylate cyclase. Several transcript variants encoding different isoforms have been found for this gene.

# **Recommended Dilutions**

**WB** 1:100 - 1:500

**IHC-P** 1:50 - 1:200

## **Immunogen Information**

**Gene ID Swiss Prot**2771
P04899

### **Immunogen**

A synthetic peptide corresponding to a sequence within amino acids 250-355 of human GNAI2 (P04899).

#### **Synonyms**

GIP; HG1C; GNAI2B; H\_LUCA15.1; H\_LUCA16.1; GNAI2

# Contact

www.abclonal.com


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

## Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,0.05% BSA,50% glycerol,pH7.3.

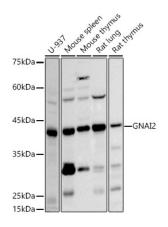

Western blot analysis of extracts of various cell lines, using GNAI2 antibody (A20930) at 1:500

dilution


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 3s.

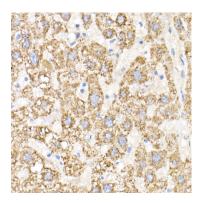

Immunohistochemistry analysis of paraffin-embedded human liver using GNAI2 Rabbit mAb (A20930) at dilution of 1:20 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

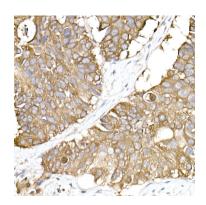

Immunohistochemistry analysis of paraffin-embedded human lung cancer using GNAI2 Rabbit mAb (A20930) at dilution of 1:20 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.