Leader in Biomolecular Solutions for Life Science

# FAM38A/PIEZO1 Rabbit pAb

Catalog No.: A4340 1 Publications



### **Basic Information**

### **Observed MW**

286kDa

### **Calculated MW**

287kDa

#### Category

Polyclonal Antibody

### **Applications**

WB,IF/ICC,ELISA

#### **Cross-Reactivity**

Human, Mouse, Rat

# **Background**

The protein encoded by this gene is a mechanically-activated ion channel that links mechanical forces to biological signals. The encoded protein contains 36 transmembrane domains and functions as a homotetramer. Defects in this gene have been associated with dehydrated hereditary stomatocytosis.

# **Recommended Dilutions**

**WB** 1:500 - 1:1000

**IF/ICC** 1:50 - 1:200

## **Immunogen Information**

**Gene ID**9780

Swiss Prot
9780

Q92508

#### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 2230-2420 of human FAM38A/PIEZO1 (NP\_001136336.2).

### **Synonyms**

ER; DHS; Mib; LMPH3; FAM38A; LMPHM6; FAM38A/PIEZO1

### **Contact**

www.abclonal.com


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

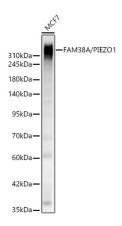

Western blot analysis of lysates from MCF7 cells, using FAM38A/PIEZO1 Rabbit pAb (A4340) at 1:1000


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25ug per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 45s.

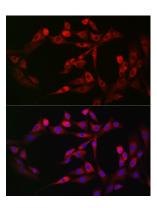


Immunofluorescence analysis of NIH-3T3 cells using FAM38A/PIEZO1 Rabbit pAb (A4340) at dilution of 1:200 (40x lens). Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

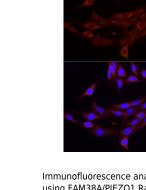

Immunofluorescence analysis of PC-12 using FAM38A/PIEZO1 Rabbit pAb (A4340) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

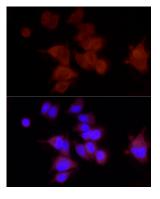

Immunofluorescence analysis of MCF-7 using FAM38A/PIEZO1 Rabbit pAb (A4340) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

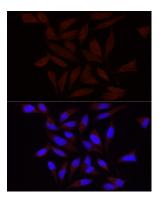

Immunofluorescence analysis of HeLa using FAM38A/PIEZO1 Rabbit pAb (A4340) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.