

A5009

Leader in Biomolecular Solutions for Life Science



GLP-2 Rabbit mAb

Catalog No.: A5009

Recombinant

2 Publications

Basic Information

Observed MW

Calculated MW

21kDa

Category

SMab Recombinant Monoclonal Antibody

Applications

IHC-P,IF/ICC,ELISA

Cross-Reactivity

Mouse,Rat

CloneNo number

ARC1143

Background

The protein encoded by this gene is actually a preproprotein that is cleaved into four distinct mature peptides. One of these, glucagon, is a pancreatic hormone that counteracts the glucose-lowering action of insulin by stimulating glycogenolysis and gluconeogenesis. Glucagon is a ligand for a specific G-protein linked receptor whose signalling pathway controls cell proliferation. Two of the other peptides are secreted from gut endocrine cells and promote nutrient absorption through distinct mechanisms. Finally, the fourth peptide is similar to glicentin, an active enteroglucagon.

Recommended Dilutions

IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200

Immunogen Information

Gene ID

2641

Swiss Prot

P01275

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 81-180 of human GLP-2 (P01275).

Synonyms

GLP1; GLP2; GRPP; GLP-1; GLP-2

Contact



www.abclonal.com

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,0.05% BSA,50% glycerol,pH7.3.

Validation Data

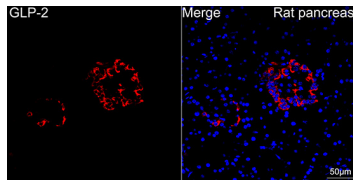

Confocal imaging of paraffin-embedded Rat pancreas tissue using GLP-2 Rabbit mAb (A5009, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform high pressure antigen retrieval with 0.01 M citrate buffer (pH 6.0) prior to IF staining.

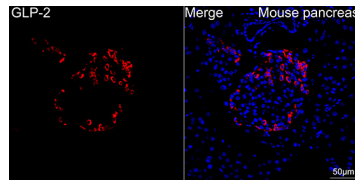

Confocal imaging of paraffin-embedded Mouse pancreas tissue using GLP-2 Rabbit mAb (A5009, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform high pressure antigen retrieval with 0.01 M citrate buffer (pH 6.0) prior to IF staining.